

3-Layer Compression Bandage System

McKESSON

McKesson 3-layer compression bandage system is a multilayer therapeutic dressing that helps deliver consistent, sustained compression for managing venous insufficiencies and other chronic venous conditions. Providing approximately 40 mmHg of resting compression, this system combines padding, a cotton-blend layer, and a cohesive compression bandage to create a secure, thin, non-slip dressing that can remain in place for up to seven days. Intended for use on intact skin and suitable for patients with an ankle brachial pressure index (ABPI) of 0.8–1.3, it supports effective long-term compression therapy with easy, hand-tear application.

Features and Benefits

Three-Layer System

Layer 1: Padding Bandage — Provides cushioning and redistributes pressure to support the management of venous insufficiencies

Layer 2: Cotton-Blend Bandage — Smooths the padding bandage to prepare the limb for even compression

Layer 3: Compression Bandage — Delivers therapeutic compression (40 mmHg) and helps keep the dressing securely in place for up to 7 days

Compression Performance

Sustained 40 mmHg Compression — Provides average resting compression of approximately 40 mmHg for effective therapy

Supports Long-Term Venous Therapy — Designed for chronic venous diseases, thromboembolic conditions, edema, and related indications

Suitable for ABPI 0.8–1.3 — Appropriate for patients within this assessment range only

Application & Usability

No Scissors Required — Bandages can be torn by hand for quick and convenient application

Air-Permeable — Allows breathability for patient comfort

Non-Slip System — When used together, layers form a secure, thin, non-slip compression dressing

Non-Sterile & Single Use — Designed for one-time application; do not reuse

Not made with natural rubber latex



3-Layer Compression Bandage System

McKESSON

Multilayer compression system (40 mmHg) to be used on intact skin

Layer 3: Compression


Provides compression (40 mmHg) and helps keep the bandage in place for up to 7 days

Layer 2: Cotton Blend

Smooths the padding bandage

Layer 1: Padding

Redistributes pressure to aid in the management of venous insufficiencies

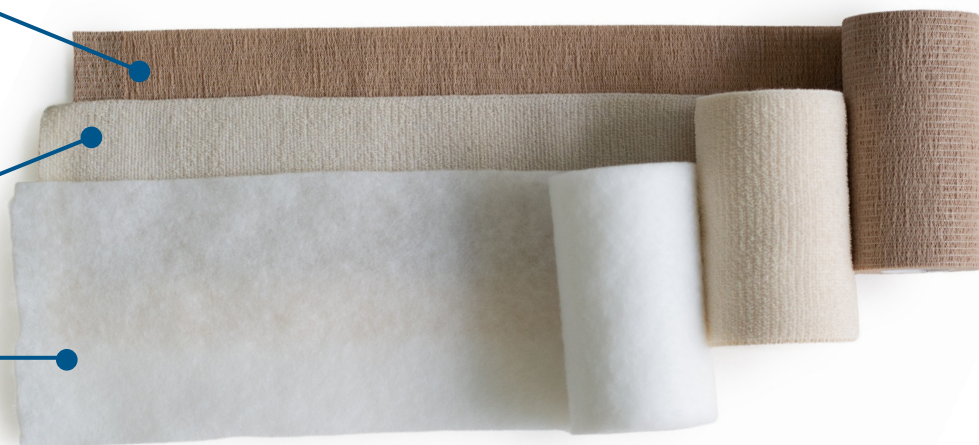

Ordering Information

E1 #	MFR #	DESCRIPTION	QTY
1278961	16-3000	3-Layer Compression Bandage System	1/kt, 8 kts/cs
		Layer 1: Padding, 4 in x 4 yds (10 cm x 3.6 m) Redistributes pressure to aid in the management of venous insufficiencies	
		Layer 2: Cotton Blend, 4 in x 6.4 yds (10 cm x 5.8 m) Smooths the padding bandage	
		Layer 3: Compression Bandage, 4 in x 6.9 yds (10 cm x 6.3 m) Provides compression (40 mmHg) and helps keep the bandages in place for up to 7 days	

Satisfaction Guaranteed

For complete details, please visit
mms.mckesson.com/mckesson-brands.

Never Compromise on Care™

With McKesson Brands, you'll never have to compromise. We deliver quality products at competitive prices, backed by our 100% satisfaction guarantee. We have decades of experience in the market, and we offer you support and work to help ensure you have the products you need, when you need them—making it even easier for you to provide the best care for your patients.

The product information contained in this document, including the product images and additional product materials, was collected from various supplier sources. All product claims and specifications are those of the product suppliers and have not been independently verified by McKesson Medical-Surgical or its affiliates ("McKesson"). McKesson is not responsible for errors or omissions in the product information.

The properties of a product may change or be inaccurate following the posting or printing of the product information in the document, either in the print or online version. Caution should be exercised when using or purchasing any products from McKesson's online or print documents by closely examining the product packaging and the labeling prior to use.

Due to product changes, information listed in this document is subject to change without notice. This information is placed solely for your convenience in ordering and McKesson disclaims all responsibility for its completeness and accuracy, whether or not the inaccuracy or incompleteness is due to fault or error by McKesson.

McKesson and the McKesson logo are registered trademarks of McKesson Corporation.

© 2026 McKesson Medical-Surgical Inc.

Rev. 00 03/26